

Lightning-Link® Rapid Biotin Conjugation Kit (Type B*)

*Optimized for assays in which conjugate is captured by a streptavidin-coated plate

Applicable to: 371-0030 3 x Ab labelings (up to 20ug)
371-0010 3 x Ab labelings (up to 200ug)

371-0005 1 x Ab labeling (up to 200µg)
371-0015 1 x Ab labeling (up to 2mg)

Release 2 © EXPEDEON. 11/09/2018

INTRODUCTION

Lightning-Link® conjugation technology works by targeting amine groups (e.g. lysines) and is widely used to label antibodies.

The Lightning-Link® Biotin kit allows Biotin conjugations to be set up in less than 30 seconds, simply by adding a solution of the antibody to a lyophilized mixture containing a proprietary activated Biotin ligand (Figure 1).

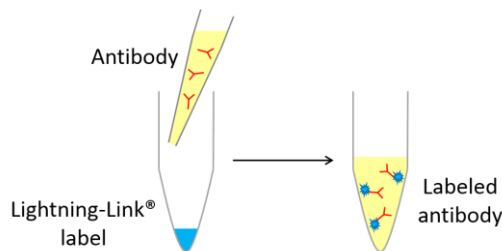


Figure 1. Lightning-Link® antibody conjugation

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, Lightning-Link® technology can be used to label both small (e.g. 10 µg) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%.

Directly labeled primary antibodies are advantageous as they eliminate the need for secondary reagents in immunoassay procedures, thus removing a tedious extra cycle of incubation and wash steps in applications, such as ELISA and Western blotting.

Lightning-Link® technology can also be used to label proteins, peptides and other biomolecules (See 'the What can I label using Lightning-Link®?' Section of the protocol for further information).

KIT CONTENTS

- 1 or 3 glass vial(s) of Lightning-Link® Rapid mix
- 1 vial of LL Rapid Modifier reagent
- 1 vial of LL Rapid Quencher reagent

SHIPPING CONDITIONS

The kit is shipped at ambient temperature in a tamper-evident polypropylene container. Store the kits at -20°C upon receipt.

Please note that the modifier and quencher FD after initial thawing can be stored at either 4°C or -20°C.

BUFFER CONSIDERATIONS

Please see the below table for recommended buffer conditions and components:

Buffer Components	
pH	6.5-8.5
Amine free buffer (e.g. MES, MOPS, HEPES, PBS)	✓
Non-buffering salts (e.g. sodium chloride)	✓
Chelating agents (e.g. EDTA)	✓
Sugars	✓
Glycerol	<50%
Thiomersal / Thimerosal	✗
Merthiolate	✗
Sodium Azide ¹	<0.1%
BSA ^{1,2}	<0.1%
Gelatin ^{1,2}	<0.1%
Tris	<50mM
Glycine	✗
Proclin	✗
Borate buffer	✓
Nucleophilic components (Primary amines e.g. amino acids or ethanolamine and thiols e.g. mercaptoethanol or DTT)	✗

¹ Please note that individually the concentrations shown should not affect the reaction. However in combination with additional compounds that are not recommended above a certain concentration, the reaction may be affected.

² If intending to use this kit for immunohistochemistry, it is recommended that there be no gelatin or BSA present.

AMOUNT AND VOLUME OF ANTIBODY

Prod. Code	Recommend amount of antibody	Maximum amount of antibody	Maximum Conjugation volume
371-0030	10µg	20µg	10µl
371-0005	100µg	200µg	100µl
371-0015	1mg	2mg	1ml

Superior conjugates are normally generated using the recommended amount of antibody. Using the maximum amount of antibody will still generate quality conjugates.

Antibodies less than 1mg/ml can still be used to generate good conjugates provided the maximum conjugation volume is not exceeded. Adding less than the recommended maximum amount of antibody may result in unbound label post conjugation. This excess label will be deactivated by the quencher and removed during the first wash step of any application. We would recommend that antibodies below 0.5mg/ml are concentrated prior to use.

Please contact our technical support team for more advice.

INSTRUCTIONS

Setting up your conjugation reaction

1. Before you add antibody to the Lightning-Link® Rapid mix, add 1 μ l of LL Rapid Modifier reagent for each 10 μ l of antibody to be labeled. Mix gently.
2. Remove the screw cap from the vial of Lightning-Link® Rapid mix and pipette the antibody sample (with added LL Rapid Modifier) directly onto the lyophilized material. Resuspend gently by withdrawing and re-dispensing the liquid once or twice using a pipette.
3. Replace cap on the vial and leave standing for 15 minutes at room temperature (20-25°C). Alternatively, and sometimes more conveniently, conjugations can be set up and left overnight at room temperature. Longer incubation times have no negative effect on the conjugation.
4. After 15min incubation, add 1 μ l of LL Rapid Quencher reagent for every 10 μ l of antibody used. The conjugate can be used after 5 minutes. The conjugates do not require purification.

STORAGE OF CONJUGATES

Your biotin conjugate can be stored at 4°C for up to 18 months. For longer storage the conjugate can be stored at -20°C with a cryoprotectant such as 50% glycerol.

The best storage conditions for any particular conjugate must be determined by experimentation.

WHAT CAN I LABEL USING LIGHTNING-LINK®?

Lightning-Link® technology works by targeting free amine groups. It can be used to label antibodies, peptides, proteins and other molecules with free amine groups.

The protocol provided here is optimized for labeling IgGs, when labeling other molecules please contact our technical support team at: www.expedeon.com/contact/

Our custom conjugation service team will also be happy to generate your conjugate for you. Please contact our sales team for a quote. www.expedeon.com/contact/

WHAT IF MY BUFFER DOESN'T FIT THE REQUIREMENTS?

The AbSelect™ purification kit range allows you to quickly and simply purify your antibody using our conjugation friendly buffers.

The appropriate kit to use depends on your particular sample. Please view our AbSelect™ range of products on our website to select the best kit.

If your antibody is already purified but its concentration is too low, you can concentrate it by using our AbSelect™ Antibody Concentration and Clean Up Kit (product code 861-0010). This kit can also be used to remove low molecular weight contaminants such as azide, Tris or glycine.

If your antibody contains BSA, you can use our AbSelect™ BSA removal kit (product code 820-0100) to purify your antibody in one simple step.

Our AbSelect™ kits are fully compatible with Lightning-Link® and are designed to work with IgG antibodies. The only exception is the Concentration and Clean Up kit (product code 861-0010) which will work with any molecule greater than 10kDa.

TECHNICAL SUPPORT

For technical enquiries get in touch with our technical support team at: www.expedeon.com/contact

Limited use license

Use and Use Restrictions: The Products are sold, and deliverables of any services are provided for the purposes of the buyer's internal in vitro research, development or educational use only, not for in vivo, or any therapeutic or diagnostic use, nor for resale, or for providing services or any other commercial use of any kind, including without limitation, for any transfer in any form (including as part of a kit) to a third party; for analysis or reverse engineering of the Product or for manufacturing. Products should only be used in accordance with any safety data sheets, guidance or protocols that we issue from time to time and are available for download from our Site. Protective clothing should be used at all times when handling our Products. Safety datasheets relating to all Products are available for download from the Site or upon request.

Expedeon grants no other license or rights under any intellectual property in respect of Products or services deliverables and in particular grants no license to use any Product or deliverables for any commercial purposes. Sale of Products or service deliverables by us or our authorised distributors are expressly conditional upon the customer's agreement with these restrictions, which the customer gives upon placing an order for Products or deliverables. If you wish to use any Product or deliverables for any purpose other than your own internal research as described above, you will require an additional licence from Expedeon. Please contact licensing@expedeon.com.